



FLORET SERIES INFANT RADIANT WARMERS

Step into Future...



Follow us on



www.bplmedicaltechnologies.com



Happier Living Everyday

INTRODUCING FLORET SERIES RADIANT WARMERS

BPL Medical introduces in India, Floret series Infant radiant warmers with a unique Skin sensor displacement alarm technology.

Now you no longer have to worry about patient overheating due to displacement of sensor probe displacement from baby skin.

Unlike the current warmers available in India, Floret series warmers incorporate a technology which detects skin sensor displaced from the baby's skin. Sensor displacement is a serious safety issue and can lead to baby overheating and fluid loss etc.

Floret series is available in several different models which will meet your operating and budget requirements.

Key Features

- 7 Segment LED displays for easy viewing and setting of baby temperature
- Three operating modes: Servo, Manual & Pre-warm.
- Audio Visual alarm: Low & high temperature, High temperature cutoff, Heater, Power & Sensor Failure, Patient alarm and Skin sensor Displacement alarm.
- 600 W Infrared quartz heater with SS parabolic reflector for faster & more accurate displacement of heat.
- APGAR timer and countdown timer display with toggle.
- 7 Watt LED lamp for observation
- Software lock for control panel.



Engineered for safety
Ergonomically design
Economical

**TECHNOLOGY MAKES
AFFORDABLE...**



DIFFERENT MODELS

FLORET 100 Series

Infant radiant warmer with non-detachable baby trolley



With
Fixed Bed & Drawer



Technical Specifications


● Modes of warmer	Three modes. Servo, Manual & Pre-warm
● Electrical specifications	80-230 VAC 2.5A, 50/60Hz Fuses 5A@220-240V, 50/60Hz
● Electronics	Microprocessor based Design
● Warmer Head	600 W Infrared Quartz Heater with SS Parabolic Reflector
● Display	Small 7 segment LED
● Temperature Sensor	Type: Thermistor Range: 26.5° to 40°C Accuracy: + / - 0.3°C Resolution: + / - 0.1°C Set Temp : 32° - 37°C
● Audio Visual Alarms	Visual and audible High Temperature Low Temperature High Temperature cut off Power failure Probe failure Check patient Heater failure Skin sensor displaced
● Observation Lamp	7 W LED lamp 6500 K
● Additional accessories	Monitor tray, SS IV Pole, Soft top, nontoxic with rigid frame
● Wheels	3 inch castor wheels with two brake wheels
● Dimensions	Length Breadth Height : 101x74x177 cms Approx. weight : 55 Kg.
● Optional Add-ons	X Ray Cassette Holder

FLORET 1000 Series

Infant radiant warmer with integrated bed and baby tray or with integrated bed and baby drawer.

Note: Floret 1000 series also available with detachable baby bassinet.

A Global Medical Technology Company

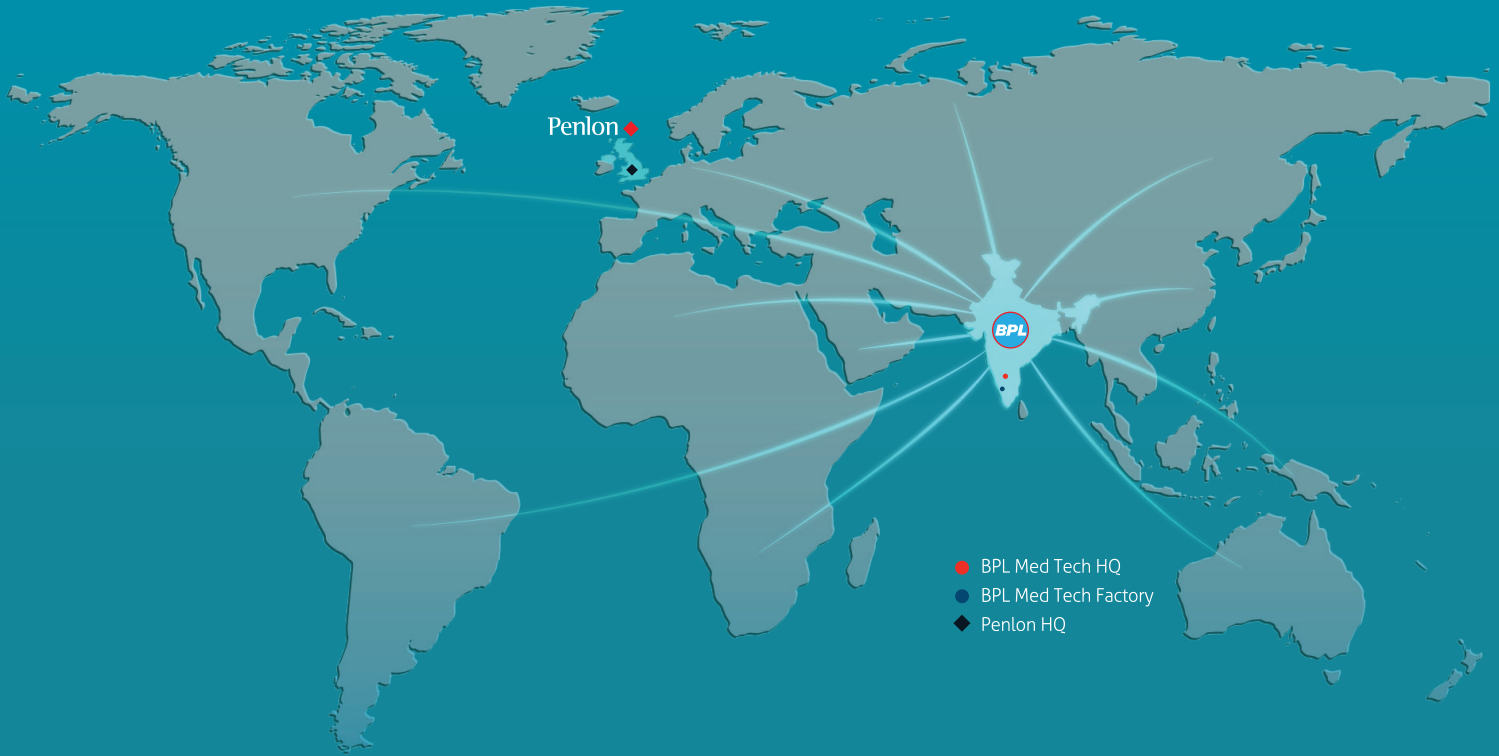

Service is Powered by:



Service Helpline
1800-425-235

CERTIFIED ISO 13485 : 2016 COMPANY

BPL MEDICAL TECHNOLOGIES PRIVATE LIMITED

Regd. Office: 11th KM, Bannerghatta Road,
Arakere, Bangalore - 560076, India.

Website: www.bplmedicaltechnologies.com

Toll Free: 1800-4252355

For Enquiries: sales.medical@bpl.in

CIN: U33110KA2012PTC067282Z



THE ECONOMIC TIMES



Happier Living Everyday